


アデノ随伴ウイルス (AAV)
Adeno-associated virus

AB220308B

1. Empty Capsid
2. Full Capsid

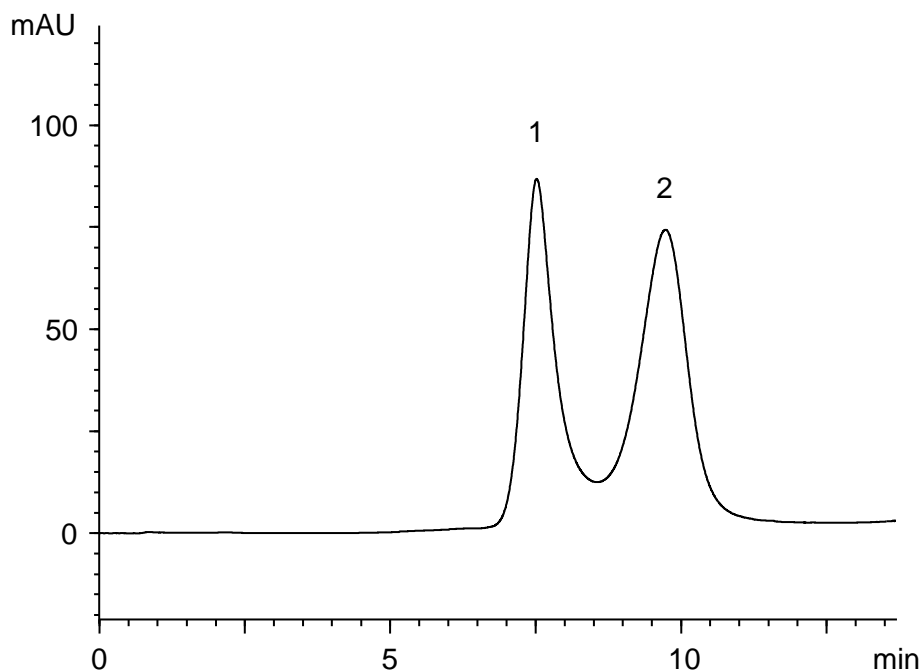

This research was supported by AMED under Grant Number JP18ae0201001.

Column	: BioPro IEX QF (5 μ m) 30 X 4.6 mm I.D.
Eluent	: A) 20 mM Bis-tris propane-HCl (pH 9.0) B) 20 mM Bis-tris propane-HCl containing 0.5 M $(\text{CH}_3)_4\text{NCl}$ (pH 9.0) 5%B (0-0.2 min), 20-45%B (0.2-10 min)
Flow rate	: 0.5 mL/min
Temperature	: 25°C
Detection	: FLS at Ex. 280 nm, Em. 348 nm
Injection	: 2 μ L
Sample	: AAV2 (2.59×10^{12} vg/mL)