




Doxepin hydrochloride

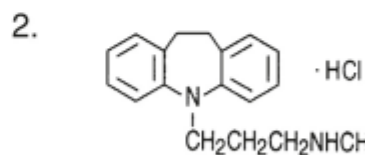

Desipramine hydrochloride

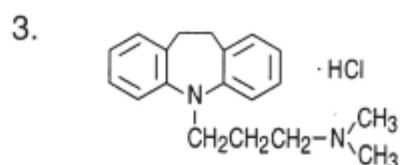

Imipramine hydrochloride

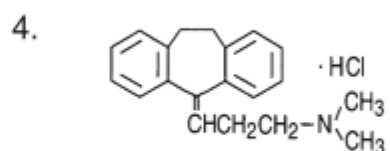

Amitriptyline hydrochloride

Column	: YMC-Pack Pro C18 (5 $\mu$ m, 12nm) 150 $\times$ 4.6 mm I.D.
Eluent	: acetonitrile/50mM NH <sub>4</sub> H <sub>2</sub> PO <sub>4</sub> (30/70)
Flow rate	: 1.0 mL/min
Temperature	: 37 °C
Detection	: UV at 215 nm, 0.13 AUFS
Injection	: 2 $\mu$ L (0.125 mg/mL)