


Application Note

YMC
EUROPE GMBH

Galantamine**

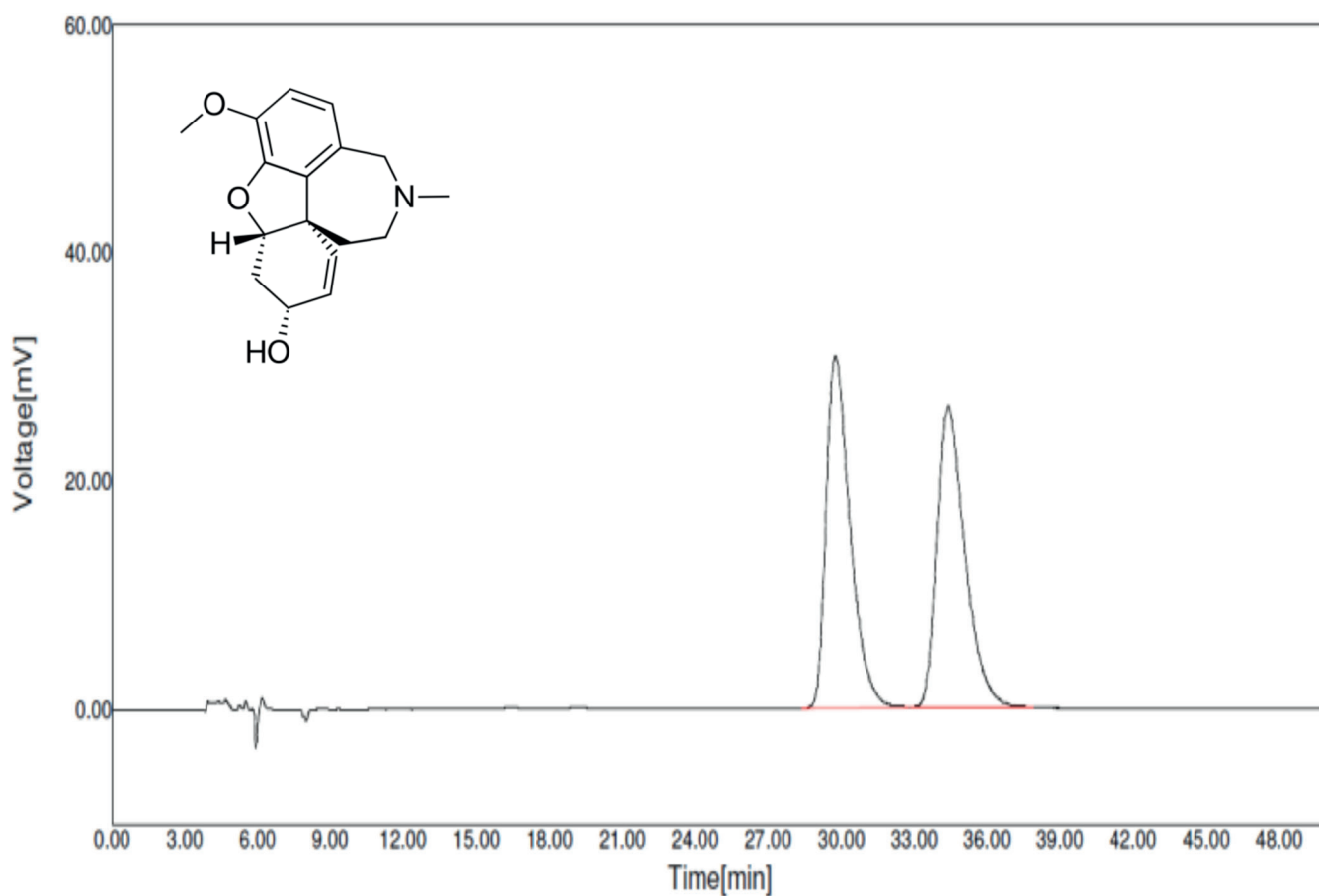

Column: CHIRAL ART Amylose-C 5 μ m (250 x 4.6 mm ID)
Part No.: KAN99S05-2546WT
Eluent: *n*-hexane/2-propanol/propionic acid/diethylamine (90/9/1/0.2)
Flow rate: 0.8 mL/min
Temperature: 15°C
Detection: UV at 289 nm
Injection: 20 μ L (0.5 mg/mL in ethanol)

**Application data by courtesy YMC India Pvt. Ltd.