


Application Note



Hexaconazole

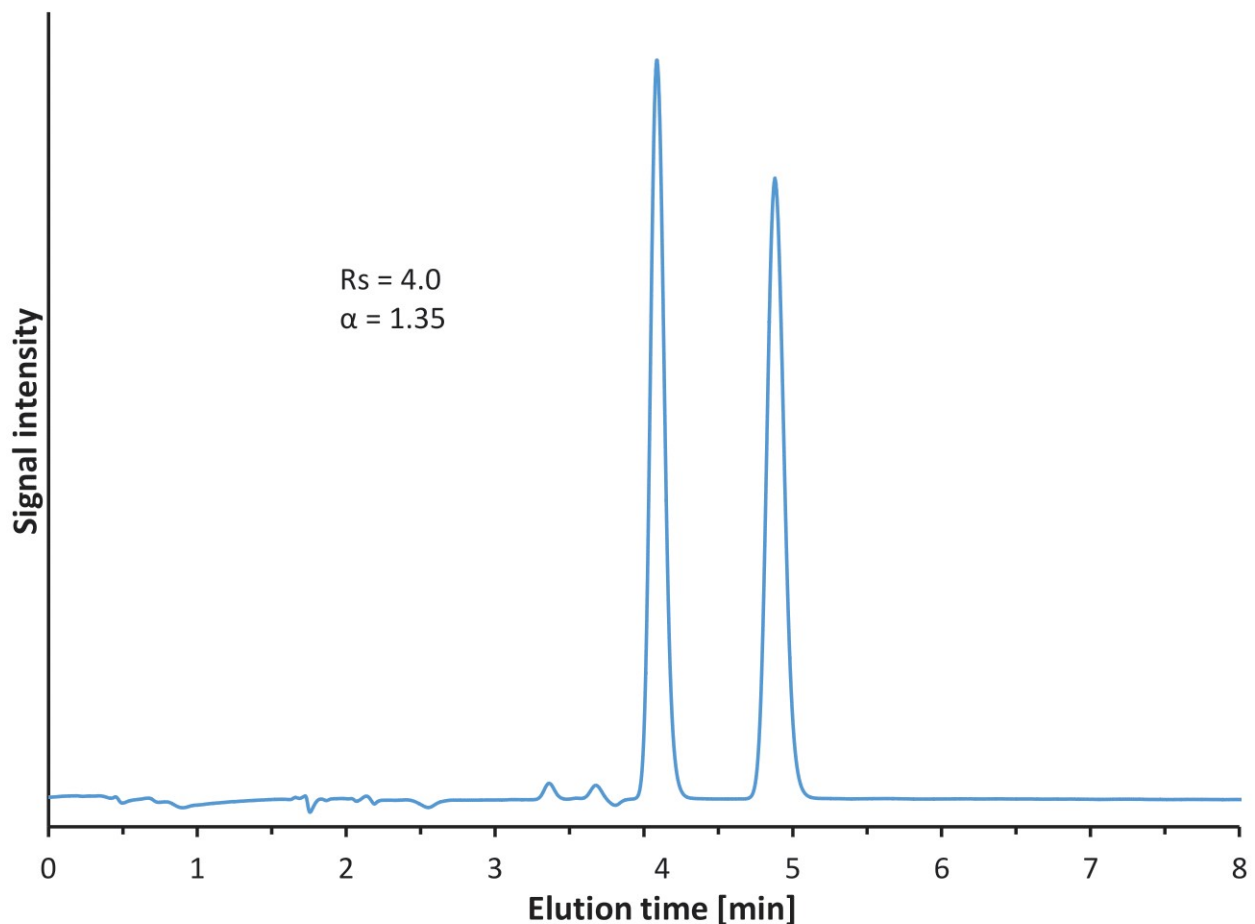

Column: CHIRAL ART Cellulose-SZ (5 μ m) 150 x 4.6 mm ID
Part No.: KSZ99S05-1546WT
Eluent: 20 mM NH_4HCO_3 / acetonitrile / diethylamine (40/60/0.1)
Flow rate: 1 mL/min
Temperature: 25 $^\circ\text{C}$
Detection: UV at 220 nm
Injection: 5 μL
Sample: Hexaconazole, 1 mg/mL 20 mM NH_4HCO_3 / acetonitrile / diethylamine (50/50/0.1)