


1.




*o*-Terphenyl

2.

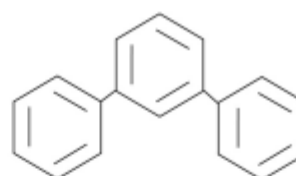

*m*-Terphenyl

3.

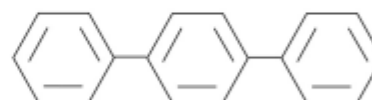

*p*-Terphenyl

Column	: YMC-Pack Pro C18 RS (5 $\mu$ m, 8nm) 150 $\times$ 4.6 mm I.D.
Eluent	: methanol/water (85/15)
Flow rate	: 1.0 mL/min
Temperature	: 37 $^{\circ}$ C
Detection	: UV at 254 nm, 0.16 AUFS
Injection	: 10 $\mu$ L (0.02 mg/mL)