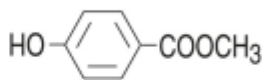
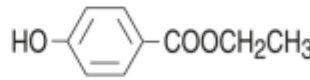
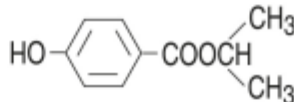
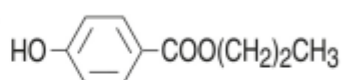
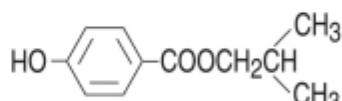
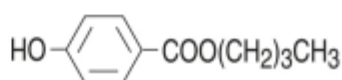


1. 
Methyl *p*-hydroxybenzoate
2. 
Ethyl *p*-hydroxybenzoate
3. 
Isopropyl *p*-hydroxybenzoate
4. 
Propyl *p*-hydroxybenzoate
5. 
Isobutyl *p*-hydroxybenzoate
6. 
Butyl *p*-hydroxybenzoate

Column	: YMC-Pack Pro C18 RS (5 μ m, 8nm) 150 \times 4.6 mm I.D.
Eluent	: methanol/water (60/40)
Flow rate	: 1.0 mL/min
Temperature	: 30 $^{\circ}$ C
Detection	: UV at 270 nm, 0.32 AUFS
Injection	: 4 μ L (0.1 mg/mL)