


1. L-Ascorbic acid (Vitamin C)
2. Orotic acid (Vitamin B13)
3. Pyrroloquinoline quinone(PQQ)
4. Pyridoxine HCl (Vitamin B6)
5. Nicotinamide

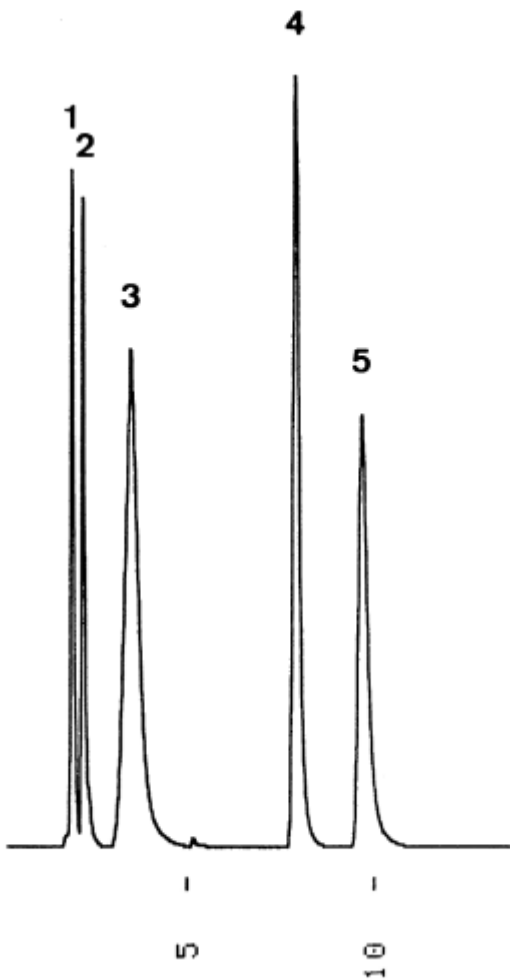

Column : YMC-Pack ODS-AM (5 μ m, 12nm)
 150 \times 4.6 mm I.D.
 Eluent : 100mM CH₃COOH/100mM CH₃COONH₄ (30/70, pH 5.1)
 Flow rate : 1.0 mL/min
 Temperature : 37 $^{\circ}$ C
 Detection : UV at 254 nm, 0.16 AUFS
 Injection : 10 μ L (0.04~0.20 mg/mL)