


- 1.p-Aminobenzoic acid (PABA)
- 2.Biotin (Vitamin H)
- 3.Riboflavin (Vitamin B₂)
- 4.Cyanocobalamin (Vitamin B₁₂)

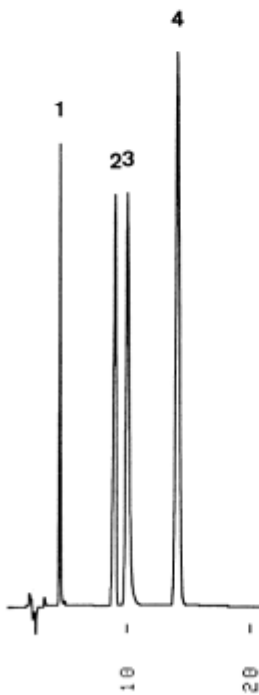

Column : YMC-Pack ODS-AM (5 μ m, 12nm)
150 \times 4.6 mm I.D.
Eluent : acetonitrile/50mM NH₄H₂PO₄ (10/90)
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 210 nm, 0.16 AUFS
Injection : 10 μ L (0.02~0.30 mg/mL)