


Application Note

YMC
EUROPE GMBH

Oxybutynin**

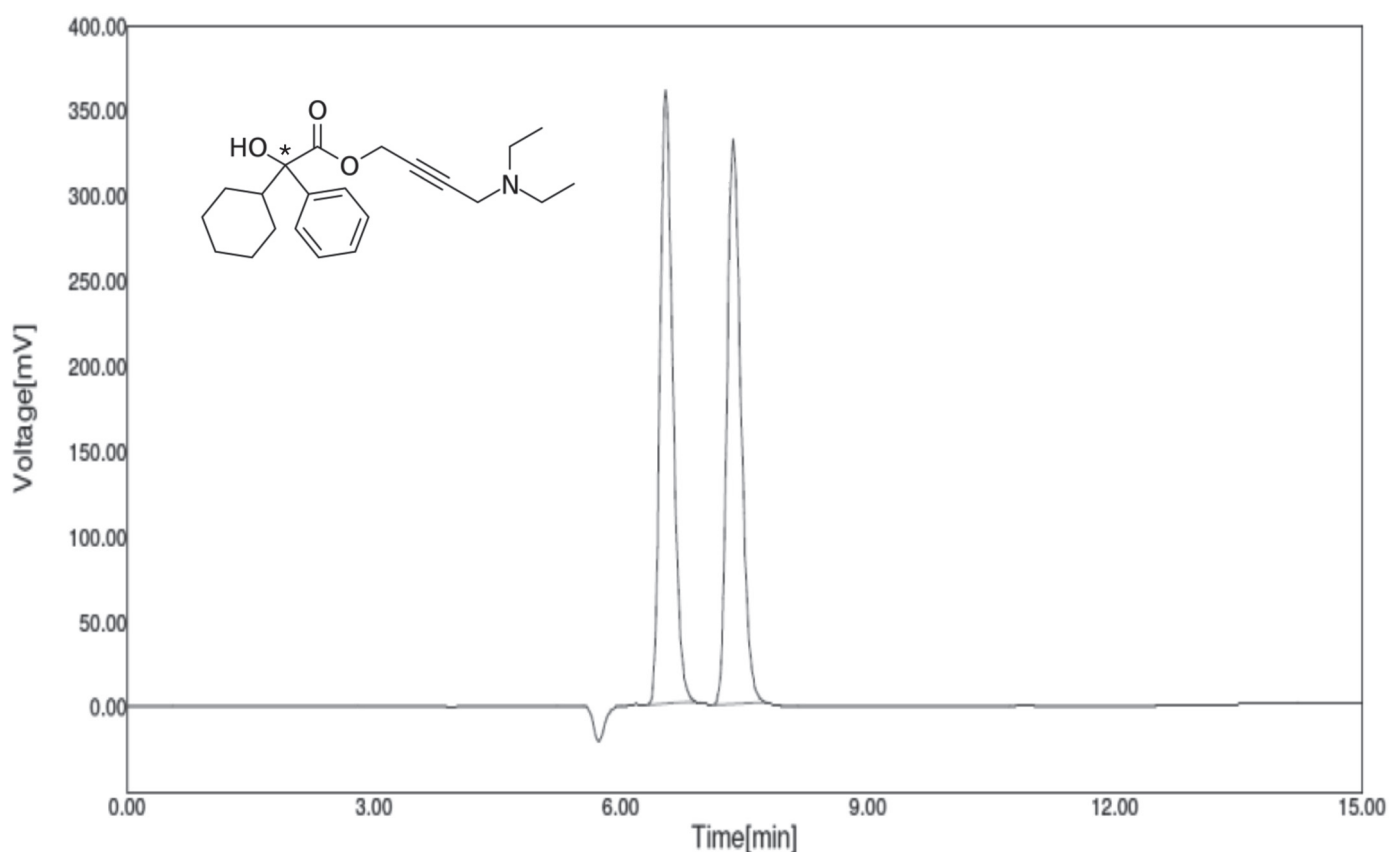

Column: CHIRAL ART Amylose-SA 5 μm (250 x 4.6 mm ID)
Part No.: KSA99S05-2546WT
Eluent: *n*-hexane/2-propanol/diethylamine (90/10/0.1)
Flow rate: 0.8 mL/min
Temperature: 25°C
Detection: UV at 220 nm
Injection: 20 μL (0.5 mg/mL)

**Application data by courtesy YMC India Pvt. Ltd.