



### Sample solution (0.4 mg/mL)

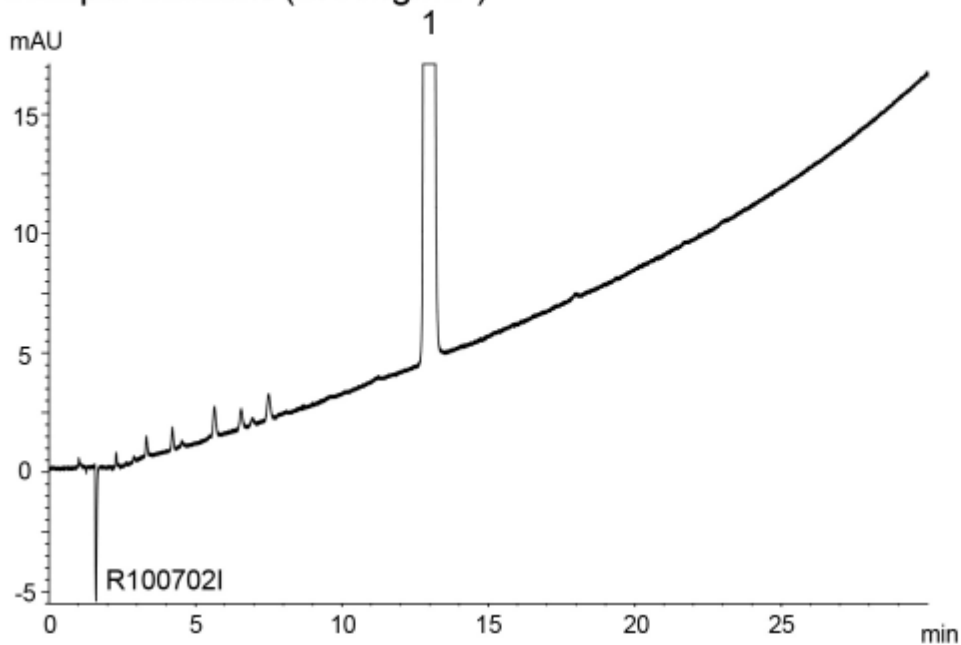


1.

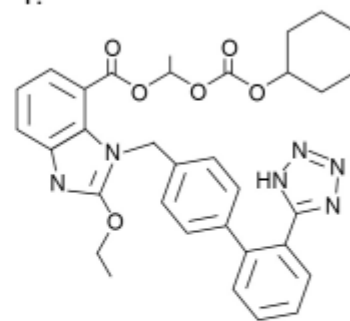

Candesartan cilexetil

### Standard solution (4 μg/mL)

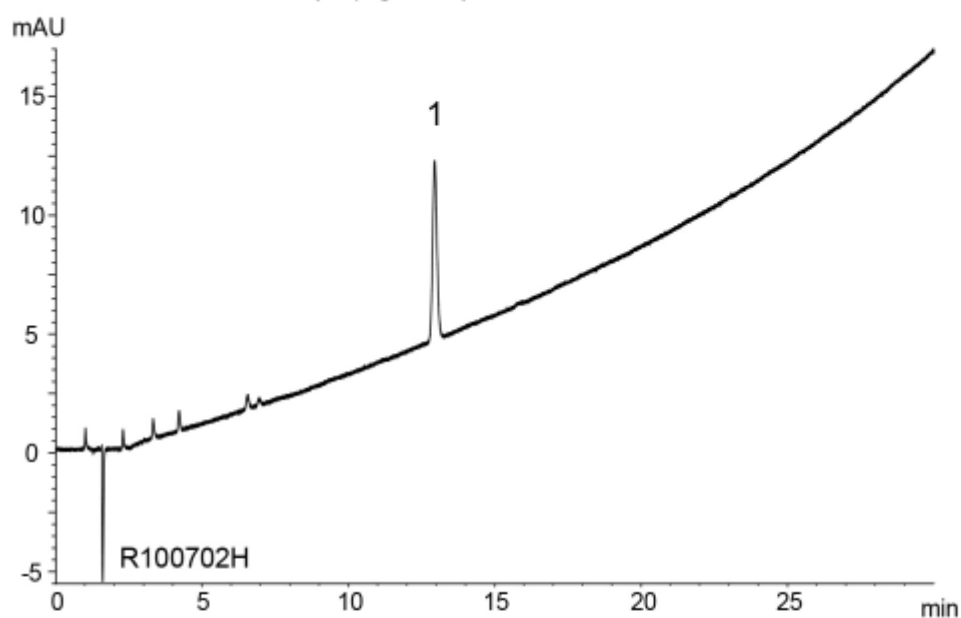

Column : J'sphere ODS-M80 (4 μm, 8 nm)  
150 X 4.0 mm I.D.  
Eluent : A) acetonitrile/water/acetic acid (57/43/1)  
B) acetonitrile/water/acetic acid (90/10/1)  
0-100%B (0-30 min)  
Flow rate : 0.8 mL/min  
Temperature : 25°C  
Detection : UV at 254 nm  
Injection : 10 μL

(日本薬局方収載原案記載条件; 純度試験)